


Simple & accurate  
**Titrating Automatic Burette**



Dynalon<sup>®</sup>  
Labware



# Available in 5 sizes from 5 to 50mL Titrating Automatic Burette

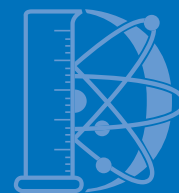

This automatic titrating burette is a universal aid in the laboratory, electroplating industry, water treatment plants, for field analyses and numerous other applications.

Automatic Titrating Burette Prominent features include:

- Self zeroing
- Measured "to deliver"
- Burette made of Schott graduated borosilicate glass with blue ring graduations and Schellbach stripe for easy to read measurements
- Graduations in accordance with Category B (Tolerances according to DIN ISO 385)
- Exact dispensing of precise discharge jet or drop via the micro screw
- Left or right hand operation for ease of use
- Sturdy stand for support
- 1 liter low density polyethylene (LDPE) chemically resistant bottle



	Burette Cap. mL	Bottle Cap. mL	Division mL	Precision mL	Height cm
896954					
-1005	5	500	0.05	±0.015	48
-1010	10	500	0.05	±0.030	57
-1015	15	500	0.10	±0.045	49
-1025	25	1000	0.10	±0.045	61
-1050	50	1000	0.10	±0.045	87

